

ApexRemote

AIRBORNE PARTICLE COUNTERS



Apex represents everything that is important in a world class Environmental Monitoring System. It is the Best In Class sensor and meets all the industry standards. The Risk Mitigation features of Self-Diagnostics and SmartBracket virtually eliminate data errors.

Apex provides real time continuous data collection at a cost effective price per point offering unmatched performance, accuracy and reliability.




Apex provides all this in an enclosure that is compatible with industry standard decontamination processes.

FEATURES

- 0.2 – 10.0 μm Particle Size Range
- 1.0 CFM or 0.1 CFM Flow Rates
- 3,000 Records Data Storage
- Self-Diagnostics
- Web Browser
- Smart Bracket
- Validation Mode
- Built in Alarm Light



LIGHTHOUSE
WORLDWIDE SOLUTIONS

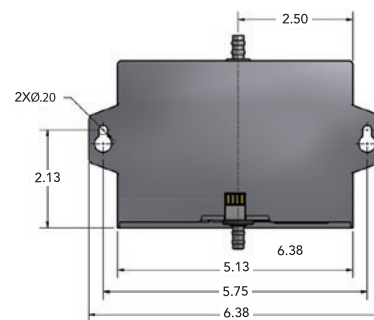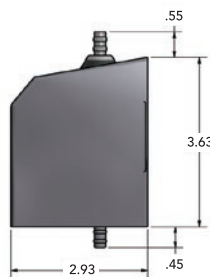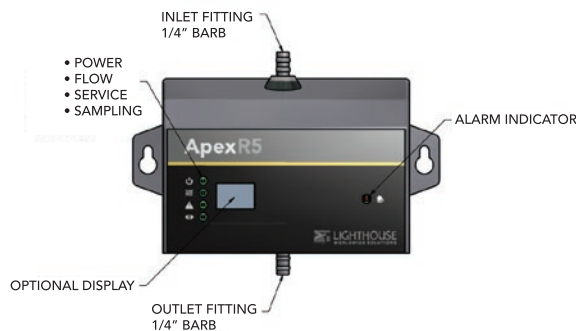

TECHNICAL DATASHEET

Apex Remote Airborne Particle Counters

Features	R3	R5	R02	R03	R05
Minimum Size Detection	0.3 µm	0.5 µm	0.2 µm	0.3 µm	0.5 µm
Standard 2 channel sizes	0.3, 0.5 µm	0.5, 5.0 µm	0.2, 0.3 µm	0.3, 0.5 µm	0.5, 5.0 µm
Optional 4 channel sizes	0.3, 0.5, 1.0, 5.0 µm	0.5, 1.0, 5.0, 10.0 µm	0.2, 0.3, 0.5, 1.0 µm	0.3, 0.5, 1.0, 5.0 µm	0.5, 1.0, 5.0, 10.0 µm
Additional Available sizes	0.3, 0.5, 0.7, 1.0, 3.0, 5.0 µm	0.5, 0.7, 1.0, 3.0, 5.0, 7.0, 10.0 µm	0.2, 0.3, 0.5, 0.7, 1.0, 2.0 µm	0.3, 0.5, 0.7, 1.0, 3.0, 5.0 µm	0.5, 0.7, 1.0, 3.0, 5.0, 7.0, 10.0 µm
Flow Rate	1 CFM (28.3 LPM)		0.1 CFM (2.83 LPM)		
Weight	1.6 lbs (0.73 kg)		1.4 lbs (0.63 kg)		
Communication	Ethernet or Serial				
Modbus Protocols	TCP or RTU or ASCII				
Data Storage Records	3,000				
Self Diagnostics	Yes				
Web Browser	Yes				
Validation Mode	Yes				
Alarm Light	Yes				
Enclosure	304 Stainless Steel				
Power	24 VDC, PoE				
External Vacuum Requirements	Nominal 1.0 CFM 18" (45.72 cm) of Hg				
Light Source	Laser Diode				
Calibration	Designed for ISO 21501-4 Compliance				
Zero Count Level	<1 count/5 minutes per ISO 21501-4				
Concentration Limits	500,000 particles/ft3 (17,150,000/m3) @ 5% coincidence loss suitable for ISO classes 1-8 cleanrooms				
Counting Efficiency	50% @ minimum channel size, 100% for particles >1.5X - 2X channel 1 particle size				
Dimensions	6.38" x 2.93" x 4.63" (16.2x7.4x11.7 cm)				

Options Available for all Models





Software LMS Express or LMS Express RT or LMS Pharma

OLED Display

SmartBracket Wall Mount Flush

SmartBracket Wall Mount Through Base

SmartBracket Wall Mount Through Wall

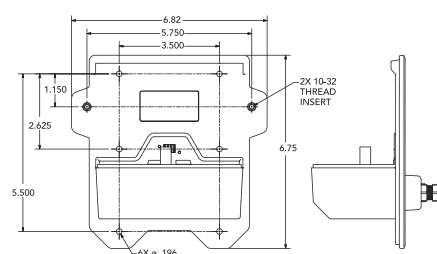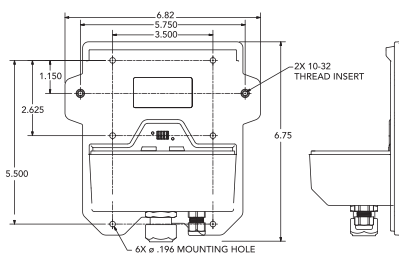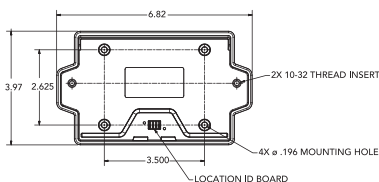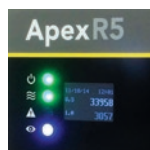

Distributed By:

laftech

Call: 1300 306 002 www.laftech.com.au

LAF Technologies Pty Ltd

Head Office: 12 Royan Place, Bayswater North, VIC 3153 Australia



Accredited Laboratory



MADE IN USA



Lighthouse Worldwide Solutions reserves the right to change specifications without notice.



LIGHTHOUSE
WORLDWIDE SOLUTIONS

www.golighthouse.com
info@golighthouse.com

1 800 945 5905 (Tollfree)
1 541 770 5905 (Outside of USA)